

Cell Marque™ Tissue Diagnostics

Claudin-4 (EP417)

Rabbit Monoclonal Primary Antibody

The claudin-family proteins are a group of tight-junction membranous proteins responsible for paracellular molecular flow between adjacent cells of the epithelium.¹ In normal kidney, claudin-4 plays a key role in paracellular ion reabsorption.² Claudin-4 is expressed in most epithelial cells yet is absent in mesothelial cells, making it of great use in characterizing epithelial malignancies.³ Immunohistochemical identification of this antigen may help identify poorly differentiated lung adenocarcinoma from lesions of mesothelial origin.⁴ Likewise, identification of claudin-4 expression by IHC may be useful in characterizing breast and endometrial lesions.^{5,6}

References:

1. Cong X, et al. Claudin-4 is required for modulation of paracellular permeability by muscarinic acetylcholine receptor in epithelial cells. *Journal of Cell Science* 2015; 128:2271-2286.
2. Hou J, et al. Claudin-4 forms paracellular chloride channel in the kidney and requires claudin-8 for tight junction localization. *Proc Natl Acad Sci USA* 2010; 107(42):18010-5.
3. Neesse A, et al. Claudin-4 as a therapeutic target in cancer. *Arch Biochem Biophys*; 2012; 524(1):64-70.
4. Ordoñez NG. Value of Claudin-4 immunostaining in the diagnosis of mesothelioma. *Anatomic Pathology* 2013; 139:611-619.
5. Abd-Elazeem MA, et al. Claudin 4 expression in triple-negative breast cancer: correlation with androgen receptors and Ki-67 expression. *Ann Diagn Pathol.* 2015; 19(1):37-42.
6. Konecny GE, et al. Claudin-3 and Claudin-4 expression in serous papillary, clear cell, and endometrioid endometrial cancer. *Gynecol Oncol.* 2008; 109(2):263-269.

Intended Use:

Claudin-4 (EP417) Rabbit Monoclonal Primary Antibody is intended for laboratory use in the detection of the claudin-4 protein in formalin-fixed, paraffin-embedded tissue stained in qualitative immunohistochemistry (IHC) testing. This product is not a stand-alone diagnostic, and cannot be used for diagnosis, treatment, prevention, or mitigation of disease.


Product Information:

Visualization: Membranous

Controls: Lung Adenocarcinoma

Dilution Range: 1:50-1:100

Associated Specialty: Pulmonary Pathology and Cytopathology

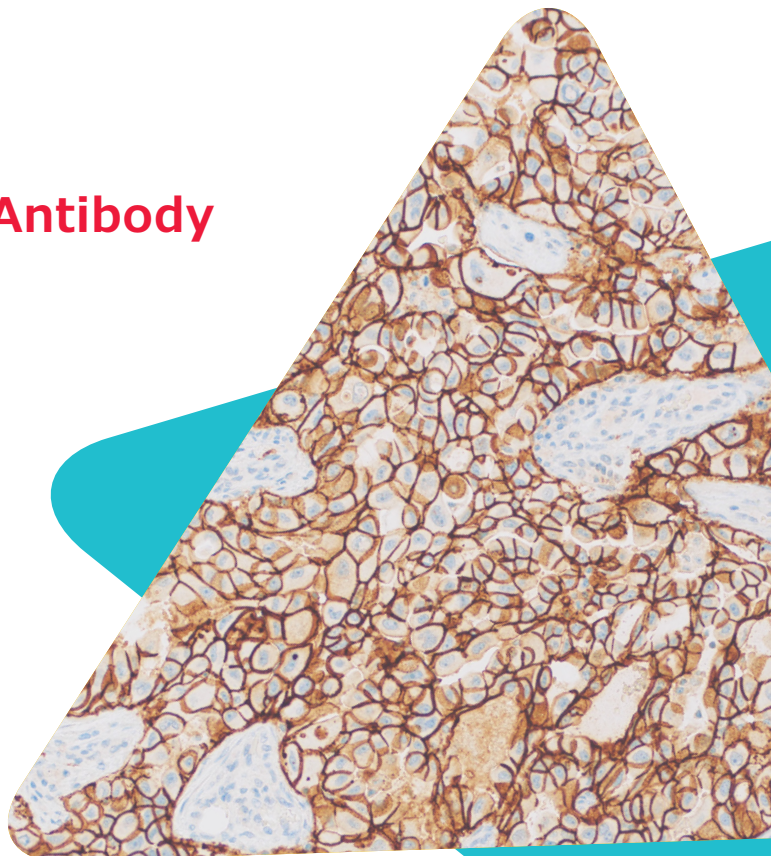


Claudin-4 on lung adenocarcinoma


Ordering Information

Description	Cat No.
0.1 mL concentrate	468R-14
0.5 mL concentrate	468R-15
1.0 mL concentrate	468R-16
1.0 mL predilute	468R-17
7.0 mL predilute	468R-18



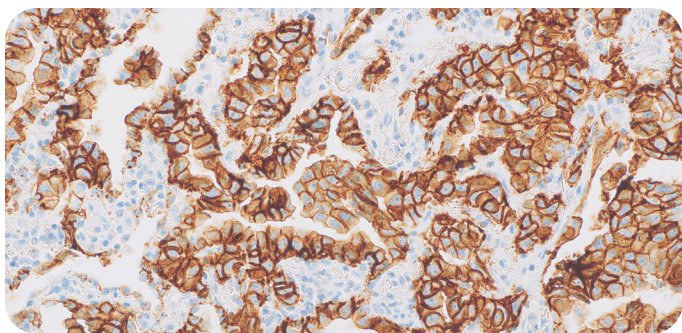

Claudin-4 on lung adenocarcinoma

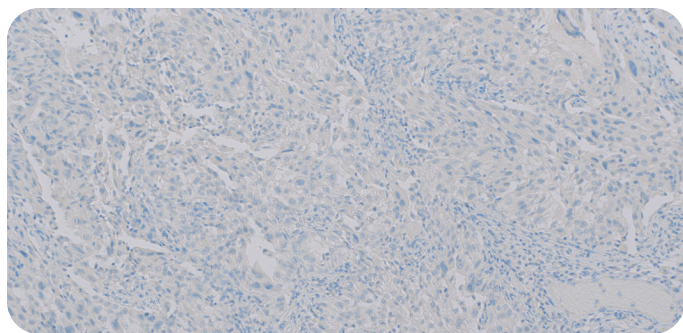

Claudin-4 on mesothelioma



USA
Toll Free: 800.665.7284
Phone: 916.746.8900
Fax: 916.746.8900
Email: service@cellmarque.com
www.cellmarque.com

CANADA
Phone: +1 916.746.8900
Fax: +1 916.746.8900
Email: international@cellmarque.com
www.cellmarque.com

MilliporeSigma
400 Summit Drive
Burlington, MA 01803

SigmaAldrich.com

© 2020 Merck KGaA, Darmstadt, Germany and/or its affiliates. All Rights Reserved. MilliporeSigma, the vibrant M, Sigma-Aldrich and Cell Marque are trademarks of Merck KGaA, Darmstadt, Germany or its affiliates. All other trademarks are the property of their respective owners. Detailed information on trademarks is available via publicly accessible resources.

MS_FL6625EN Ver. 1.0 30315 09/2020

